

Auto Refracto-Keratometer User's Manual



Please be sure to read this manual carefully before using the instrument and keep it handy for ready reference.

The Auto Refractometer is one high precision instrument of objective measuring the patient's eyes with unique optical system inside and accurate imaging analyzing and processing in Hartman technology. It's mainly used to measure the patient's diopter, including sphere power, cylinder power, optical axis, pupil distance and corneal curvature, to provide reference data for eyes' treating and eyeglasses choice. The measurement result can be displayed on screen or printed out on paper, and can also be transferred to auto phoropter (fit to RS232 interface)). If the auto phoropter can output data, this Auto Refractometer can directly print out the measured optometry data by auto phoropter.

SPECIFICATIONS

Differences	Auto Refractometer-Keratometer	Auto Refractometer
Production	Edit and embed software of refractometer and keratometer, use objective standard tool, corneal curvature standard device to adjust and calibrate	Edit and embed software of refractometer only, use objective standard device to adjust and calibrate
mode	REF 、 KER 、 and REF\KER three measurement modes optional	Default REF measurement mode
Imaging principle design	A standard ring is projected onto the patient's cornea with a group of 760nm infrared light, the corneal astigmatism and axis are calculated according to the different curvature parameters in each direction of the corneal reflex ring. Use 850nm infrared light to project the Hartmann lenslet array into the fundus of the patient's eyes, and separately calculate the parameters such as sphere, cylinder and axis according to the fundus reflection array.	Use 850nm infrared light to project the Hartmann lenslet array into the fundus of the patient's eyes, and separately calculate the parameters such as sphere , cylinder and axis according to the fundus reflection array.
Purpose	The instrument can be used to measure the curvature of the meridians of the center area about 3mm on the anterior corneal surface, that is the radius of curvature and curvature, so as to determine whether or not the cornea has astigmatism, astigmatism power and axial direction. It can also be used to measure the power of the patient, such as sphere, cylinder and axis.	used to measure the power of the patient, such as sphere, cylinder and axis.

Measurement mode	Manual measurement and automatic measurement. Aligning the measurement target at the corneal vertex, manual or automatic focus to clear state, the instrument will automatically measure the corneal curvature power. Use the same measuring method to automatically measure the powers of sphere, cylinder and axis.	Manual measurement and automatic measurement. Aligning the measurement target at the corneal vertex, manual or automatic focus to clear state, the instrument will automatically measure the powers of sphere, cylinder and axis.
Patient	Non-contact measurement, the patient opens his eyes and looks right ahead, and move the target to complete the measurement(note the replacement of the chinrest paper and the forehead rubber cleaning)	Non-contact measurement, the patient opens his eyes and looks right ahead, and move the target to complete the measurement(note the replacement of the chinrest paper and the forehead rubber cleaning)

The serial Auto Refractometer consists of optical system, mechanical transmission system, COMS image sensing system, microcomputer control system, and printer etc., are the professional instruments of providing reference data for glasses and eye diagnosis and treatment. According to measurement function, it's divided into refractive parameters measurement and refractive parameters measurement/corneal curvature measurement. According to measurement mode, it's divided into manual measurement and automatic measurement.

DISCLAIMER

1. This manual has been carefully checked to insure the contents' accuracy and perfect during compiling, however, for possible errors or omissions contained herein.
2. The company reserves the right to make changes to this product or the specifications at any time without prior notice.
3. The company owns the final interpretation to this manual.

CONTENTS

1. ELECTROMAGNETIC COMPATIBILITY GUIDE AND MANUFACTURER STATEMENT.....	1
2. SAFETY PRECAUTIONS.....	7
3. UNPACKING AND INSTALLATION	11
4. FUNCTIONS OF THE MAJOR COMPONENTS.....	13
5. MAIN TECHNICAL INDEXES.....	15
6. ENVIRONMENT TERMS.....	16
7. LCD SCREEN DISPLAY.....	17
8. MENU.....	18
9. MEASUREMENT	26
10. COMMON TROUBLE SHOOTING	31
11. PACKAGING, TRANSPORTATION, STORAGE.....	32
12. ENVIRONMENTAL PROTECTION.....	32
13. ACCESSORIES	32

1. ELECTROMAGNETIC COMPATIBILITY GUIDE AND MANUFACTURER STATEMENT

This product is in compliance with the electromagnetic compatibility regulations in this manual. To ensure compliance with these regulations, the user needs to install and use the information provided in this manual. Such as the use of non-manufacturers to provide the cable may cause the increase or decrease in the immunity of the product launch.

Warning!

1. The use of non-manufacturer supplied cables may cause an increase in the electromagnetic radiation of this product or decrease the immunity.
2. Portable or mobile radio frequency communication equipment should not be used closer to any part of the auto refractometer than the recommended isolation by distance, including the cable.
3. In addition to the transducer and cable in sale as spare parts of components from the original equipment or system manufacturer, the use of other accessories, transducers and cables may cause an increase in the device or system to launch or decrease in immunity.
4. The device or systems should not be close to or stacked up with other devices, and if you have to approach or stack, it should be observed to verify the normal operation of its use.
5. The other accessories, transducer or cable to be used together with the device and system, it may cause an increase in the device or system to launch or decrease in immunity.

1.1 Electromagnetic Emission Guide and Manufacturer Statement(Form1)


Guide and manufacturer's statement—Electromagnetic emission		
[Prototype ARK-B] expected to be used in the electromagnetic environment of the following requirements, buyers and users should ensure that it is used in this electromagnetic environment		
Launching Test	Conformity	Electromagnetic Environment—Guide
Radio frequency emission GB 4824	Group 1	[Prototype ARK-B] Radio frequency energy to be used for internal function only. Therefore, its RF emission very low, and the possibility of interference in the electronic device is very small.
Radio frequency emission GB 4824	Class B	[Prototype ARK-B] Applicable for all of the facilities in use, including the home and the direct connection of residential public low voltage power supply network.
Harmonic emission GB 17625.1	Not applicable	
Voltage fluctuation/Flicker emission GB 17625.2	Not applicable	

1.2 Electromagnetic Immunity Guide and Manufacturer Statement (Form2)

Guide and manufacturer's statement—Electromagnetic immunity			
[Prototype ARK-B] expected to be used in the electromagnetic environment of the following requirements, buyers and users should ensure that it is used in this electromagnetic environment			
Immunity Test	IEC60601 Test Level	Meet Level	Electromagnetic Environment—Guide
Electrostatic discharge GB/T 17626.2	±6kV contact discharge ±8kV air discharge	±6kV contact discharge ±8kV air discharge	The ground should be wood, concrete or ceramic tile, if the ground is covered with synthetic material, the relative humidity should be at least 30%

Electric fast transient pulse group GB/T 17626.4	±2kV power line ±1kV input/output line	±2kV power line	Network power supply should have a typical commercial or hospital environment in the use of quality
Surge GB/T 17626.5	±1kV line to line ±2kV line to ground	±1kV line to line ±2kV line to ground	Network power supply should have a typical commercial or hospital environment in the use of quality
Power input line voltage dips, short interruptions and voltage variations GB/T 17626.11	<5%U _t , last 0.5 cycle (Above U _t , >95% sag) 40% U _t , last 5 cycle (Above U _t , 60% sag) 70% U _t , last 25 cycle (Above U _t , 30% sag) <5% U _t , last 5s (Above U _t , >95% sag)	<5%U _t , last 0.5 cycle (Above U _t , >95% sag) 40% U _t , last 5 cycle (Above U _t , 60% sag) 70% U _t , last 25 cycle (Above U _t , 30% sag) <5% U _t , last 5s (Above U _t , >95% sag)	Network power supply should have a typical commercial or hospital environment in the use of quality. If the users need [Prototype ARK-B] to continuously run during power supply interruption, then it's recommended the [Prototype ARK-B] is powered by a constant power supply or battery
Power frequency magnetic field (50Hz) GB/T 17626.8	3A/m	3A/m	The power frequency magnetic field should have the characteristics of the power frequency magnetic level in a typical commercial or hospital environment
Note: U _t refers to the AC network voltage before applying the test voltage.			

1.3 Electromagnetic Immunity Guide and Manufacturer Statement (Form3)

Guide and manufacturer's statement—Electromagnetic immunity			
<p>[Prototype ARK-B] expected to be used in the electromagnetic environment of the following requirements, buyers and users should ensure that it is used in this electromagnetic environment</p>			
Immunity Test	IEC60601 Test Level	Meet Level	Electromagnetic Environment—Guide
<p>Radio frequency transmission</p> <p>GB/T 17626.6</p> <p>Radio frequency radiation</p> <p>GB/T 17626.3</p>	<p>3 V (effective value)</p> <p>150 kHz ~ 80 MHz</p> <p>3 V/m</p> <p>80 MHz ~ 2.5 GHz</p>	<p>3V (effective value)</p> <p>3 V/m</p>	<p>Portable or mobile radio frequency communication equipment should not be used closer to any part of [Prototype ARK-B] refractometer than the recommended isolation by distance, including the cable. The distance should be calculated with the corresponding formula of the transmitter frequency. The recommended isolation distance:</p> <p>$d=1.2\sqrt{P}$ </p> <p>$d=1.2\sqrt{P}$ 80MHz~800MHz</p> <p>$d=2.3\sqrt{P}$ 800MHz~2.5GHz</p> <p>In formula: P — Maximum output rated power of the transmitter provided by the manufacturer, unit for Watt(W)</p> <p>d—Recommended isolation distance, unit for meter(m).</p> <p>The electric field intensity of fixed radio frequency transmitter is determined by the investigation ^a of electromagnetic field, in each frequency range ^b should be lower than Meet Level. Interference may occur near the devices marked with the following items.</p>

Note 1: at 80MHz and 800MHz frequency point, use the formula for higher frequency bands

Note 2: these guidelines may not be suitable for all cases, because the electromagnetic propagation is influenced by the absorption and reflection of buildings, objects and human bodies.

^a fixed transmitter, such as wireless(cellular/cordless) telephone and ground mobile radio base station, amateur radio, Am and FM radio and television broadcasting, etc. the electric field intensity can not be accurately predicted in theory. In order to evaluate the electromagnetic environment of a fixed RF transmitter, the survey of electromagnetic field should be considered. If the electric field intensity measured where [**Prototype ARK-B**] place is higher than above applicable RF Meet Level, [**Prototype ARK-B**] should be observed to verify whether it can work normally. If abnormal performances happen, the supplementary measures may be necessary, such as re-adjust the direction or position of[**Prototype ARK-B**]

^b in the entire frequency range of 150 kHz~80 MHz, the electric field intensity should be less than 3V/m.

1.4 The Recommended Isolation Distance Between Portable and Mobile Radio Frequency Communication Equipments and [**Prototype ARK-B**] (Form 4)

The recommended isolation distance between portable and mobile radio frequency communication equipments and [**Prototype ARK-B**]

[**Prototype ARK-B**] expected to be used in the electromagnetic environment of the radio frequency radiation disturbance controlled. According to the maximum output rated power of communication equipment, the buyer or user may prevent the electromagnetic interference by maintaining a minimum distance to be recommended as following items between the portable and mobile radio frequency communication equipment(transmitter) and [**Prototype ARK-B**]

Maximum output rated power of transmitter: W	Isolation distance of different frequency of transmitter/m		
		150kHz~ 80MHz	80MHz~ 800MHz


	$d=1.2\sqrt{P}$	$d=1.2\sqrt{P}$	$d=2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

To the maximum output rated power of transmitter that not listed in the above forms, **d** is recommended as isolation distance, unit for meter(m), the formula in the frequency column of the corresponding transmitter is available, here **p** is the maximum output rated power of transmitter provided by the manufacturer, unit for Watt(W).

Note 1: at 80MHz and 800MHz frequency point, use the formula for higher frequency bands

Note 2: these guidelines may not be suitable for all cases, because the electromagnetic propagation is influenced by the absorption and reflection of buildings, objects and human bodies.

2. SAFETY PRECAUTIONS

Safety Signs and Instructions  Be Careful, remove the screw before installation

2.1 Operation

- 2.1.1 Don't optionally open and touch the inside parts of the instrument, it may cause an electric shock or the system may malfunction.
- 2.1.2 Please keep this instrument ground connection well to avoid possible injury to people or the instrument damaged.
- 2.1.3 Don't touch the screen(resistance touch-screen, contact and hold it for about 0.1 second) in too strong strength, it may damage the screen.
- 2.1.4 Don't put the instrument at the place of direct sunlight or too strong illumination, it may affect the measuring precision. It's strongly suggested to be used indoors or in darkroom.
- 2.1.5 Don't use the instrument in a hot, humid or dusty environment. Such environments cause bad influences to the instrument.
- 2.1.6 If you want to connect this instrument to other instrument, please follow our local agent's instructions.
- 2.1.7 In cold room, when temperature suddenly rise, dew maybe appear on the protection glass of measuring window or internal optical parts. In case this happen, it can be used till the dew disappear.
- 2.1.8 Keep the measuring window lens clean at all time. The dust and other substances may cause error in measuring or affect the measuring precision.
- 2.1.9 If you encounter any abnormal conditions, such as smoking or strange smells, turn off the instrument and pull out the power cord immediately. Contact the local experts/agent or original manufacturer to check and repair, you can use till the trouble is absolutely removed.
- 2.1.10 The use of materials that directly contact with the skin part: During operating the instrument, it should be used to separate the instrument from the patient's touch part with medical non-woven fabric(size for 8cm x 8cm), to avoid the direct contact with the surface of the instrument.

2.2 In Storage

- 2.2.1 Don't store the instrument in a place where it may get wet or where poisonous gas or liquid is stored.
- 2.2.2 Be sure to store the instrument in a place away from direct sunlight and with the specified temperature and humidity.

2.3 In Transference

- 2.3.1 During carrying the instrument, please take great care to avoid colliding and falling. Sudden or strong impact may damage the instrument or performances.

- 2.3.2 Before carrying, please turn off the machine and lock tightly the sliding body. During carrying, please catch the bottom tightly by two hands.

2.4 After Using












- 2.4.1 If the instrument won't be used for a long time, disconnect the power cable from the wall-outlet. It may cause a fire.
- 2.4.2 When the instrument is not used, turn the power off and put the dust cover on. Keeping the machine in electricity supplying will reduce the use life of the instrument. If the instrument is not covered for a long time, dust may affect the measuring accuracy.



2.5 In Maintenance

- 2.5.1 It's one high precision optical instrument and need to be calibrated regularly.
- 2.5.2 Please fill the lube to the sliding parts regularly at the experts' guide.
- 2.5.3 Be sure to replace the fuse after disconnecting the power cord from the power inlet and use the specified fuse for replacement. Otherwise, it may cause a fire.
- 2.5.4 In case the instrument breaks down, it must be checked and repaired by the specified experts who know this instrument very well, or contacting the local authorized agent or original manufacturer. Open and repair the instrument by oneself, the agent or manufacturer don't be in charge of the consequence.
- 2.5.5 This instrument for non sterile medical device. Daily cleaning and disinfection of device by end user. Please use a soft cloth or sponge, wet cloth or detergent to clean the device. Don't use alcohol, water, benzene and other organic compounds to clean the surface of the instrument, to avoid damage to the device. The measuring window is often cleaned by a soft cloth to remove dust to maintain the accuracy of the measurement.
- 2.5.6 Determination of disinfection method:
According to the requirements of "environmental and object surface disinfection" in the hospital disinfection and sanitation standard of GB 15982-2012, the general components of the device are cleaned in time, the chin-rest and forehead rest and operation lever are demanded to be disinfected in middle level disinfection.
In accordance with the regulations of the WS/T 367-2012 medical institutions disinfection technical specifications, the chin-rest and forehead rest and operation lever are demanded to be disinfected by using alcohol disinfectant, and use 75% (volume ratio) ethanol solution to wipe the surface of the object.

2.6 Safety Symbol

The International Electrotechnical Commission(IEC) has established a set of symbols for medical electronic equipment which classify a connection or warn of any potential hazards. The classifications and symbols are shown below.

	<p>I and O on power switch represent ON and OFF respectively</p>
	<p>Type B Isolated patient connection</p>
	<p>This symbol identifies a safety note. Ensure you understand the function of this control before using it. Control function is described in the appropriate User's or Service Manual.</p>
	<p>It indicates the year of manufacture and the manufacturer.</p>
	<p>Manufacturer</p>
	<p>Authorised Representative in the European Community</p>
	<p>Identifies the point where the system safety ground is fastened to the chassis. Protective earth connected to conductive parts of Class I equipment for safety purposes.</p>
	<p>Temperature Limitation</p>
	<p>Keep DRY</p>
 <p>UL60601-1 CAN/CSA C22.2 NO.601.1</p>	 <p>MEDICAL EQUIPMENT WITH RESPECT TO ELECTRIC SHOCK FIRE AND MECHANICAL HAZARDS ONLY</p>

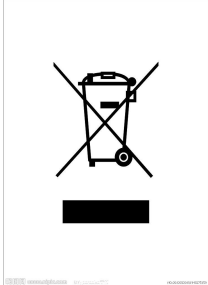	<p>ACCORDANCE WITH UL 60601-1, AND CAN/CSA C22.2 NO.601.1</p>
	<p>Disposal of your old appliance When this crossed-out wheeled bin symbol is attached to a product it means the product is covered by the European Directive 2002/96/EC. All electrical and electronic product should be disposed of separately from the municipal waste stream via designated collection facilities appointed by the government or local authorities. The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.</p>
	<p>Alternating Current</p>

3. UNPACKING AND INSTALLATION

3.1 Notices and Procedures of Taking out the Instrument

Catching the bottom and chinrest frame separately by two hands, don't catch the screen or operation lever (Figure 1)



(Figure 1)

3.2 Power Line Inlet and RS232 Interface (Figure2、 3)

Connect the spare power line with the power supply socket (RS232\USB interface connection demanded, please contact the original manufacturer or local agent)



Figure2



Figure3

Power Supply Socket: Power Switch/Socket: AC power on/off with indicator light / AC power inlet (fuse F5AL 250V inside)

Data Interface:Data Interface: USB\RS232 interface to be connected with other equipments (This refractometer and the automatic phoroptor serials connected should comply with the related electrical requirements in IEC60601-1-1-2000)

3.3 Chinrest Paper Installation

Use the specified chinrest paper (Figure 4)

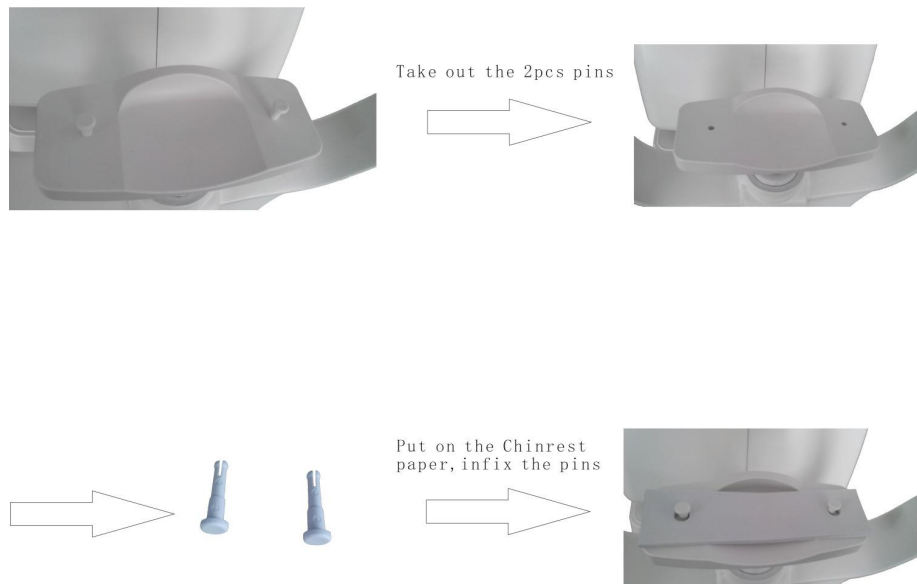


Figure 4

4. FUNCTIONS OF THE MAJOR COMPONENTS

Front (Figure 5)



(Figure 5)

LCD Screen: Monitor for measurement display

Height Adjustment Mark: The eyes' height position of the patients

Printer Cover: Press the cover to open or close

Measure Button: Performing the measurement by pressing after focusing

Joystick: Adjust the focus by moving it left/right, up/down, forward/backward

Stage Fixing Lever: Lock the sliding body

Rubber Feet: Support and adjust the instrument horizontal

Chinrest Up/Down Button: Adjust the height of the chinrest

Back (Figure 6)



(Figure 6)

Chinrest: The platform for placing the patients' chin

Forehead Rest: The place against the patients' forehead

Measuring Window: Imaging on the retina of the patients' eyes

5. MAIN TECHNICAL INDEXES

5.1 Measurement Performance Parameters

- 5.1.1 Corneal Vertical Distance(VD): 0.00mm、12mm、13.75mm、15mm
- 5.1.2 SPH: $-20.00\text{ m}^{-1}\sim +20.00\text{ m}^{-1}$ (VD=12mm, 0.01 m^{-1} 、 0.06 m^{-1} 、 0.12 m^{-1} 、 0.25 m^{-1} unit), deep myopia measurement available
- 5.1.3 CYL: $0.00\text{ m}^{-1}\sim \pm 6.00\text{ m}^{-1}$ (0.25 m^{-1} unit)
- 5.1.4 Cylinder Form: $-$ 、 $+$ 、 \pm
- 5.1.5 Axis(AX): $1^{\circ}\sim 180^{\circ}$ (1° unit)
- 5.1.6 Pupil Distance(PD): 10mm \sim 85mm(0.1mm unit)
- 5.1.7 Radius of Corneal Curvature: 5.0 \sim 10.0mm (0.01mm unit)
- 5.1.8 Corneal Power: $33.00\text{ m}^{-1}\sim 67.00\text{ m}^{-1}$ (in case that the corneal equivalent refractive power is 1.3375)
- 5.1.9 Corneal Astigmatism: $0.00\text{ m}^{-1}\sim 15.00\text{ m}^{-1}$ ($0.06\text{ m}^{-1}/0.12\text{ m}^{-1}/0.25\text{ m}^{-1}$ unit)

5.2 Other Performance Parameters

- 5.2.1 7" TFT touch screen (angle adjustable)
- 5.2.2 Printer: 57mm thermal printer
- 5.2.3 Measuring Light Energy: $<30\mu\text{W}$ (prevent injury to eyes during measuring)
- 5.2.4 Measuring Time: $<0.5\text{s}$
- 5.2.5 Minimum 2.0mm pupil can be measured. The application of cloud and mist chart technology allows the patients' eyes to look at the internal targets in a natural and comfortable situation and make the measurement more accurate
- 5.2.6 Auto tracking, auto focusing and auto measuring of end of measurement (partial model)
- 5.2.7 Electrical Power: AC100 \sim 240V, 50/60Hz
- 5.2.8 Consumption: 60AV
- 5.2.9 N.W.: 24kgs
- 5.2.10 G.W.: 26.5kgs
- 5.2.11 Dimensions: L528mm \times W260mm \times H448mm
- 5.2.12 The refractometer service life for 10 years, to ensure accurate measurement, please make the metrological verification every year

5.3 Protection Level

- 5.3.1 Product Grade: Medical apparatus and instruments grade II
- 5.3.2 Electric Shock: Level I (ground)
- 5.3.3 Electric Shock Protection Class: Class B


5.4 Device Type

- 5.4.1 Anti Electric Shock Type: Class I
- 5.4.2 Anti Electric Shock Degree: Applicable type B
- 5.4.3 Non AP device, non APG device
- 5.4.4 Running Mode: Continuous duty

6. ENVIRONMENT TERMS

- 6.1 Temperature: 10°C ~ 30°C
- 6.2 Relative Humidity: (30~75) %RH
- 6.3 Atmospheric Pressure: 86kPa ~ 106kPa
- 6.4 Altitude: <2000m
- 6.5 No strong vibration and corrosive gas around
- 6.6 No strong electromagnetic interference around
- 6.7 Brightness: <150Lx
- 6.8 The device should be placed at the specified instrument table that can rise and fall vertically
- 6.9 The device can't be used in the environment of flammable and anesthetic gas

7. LCD SCREEN DISPLAY (Figure 7)

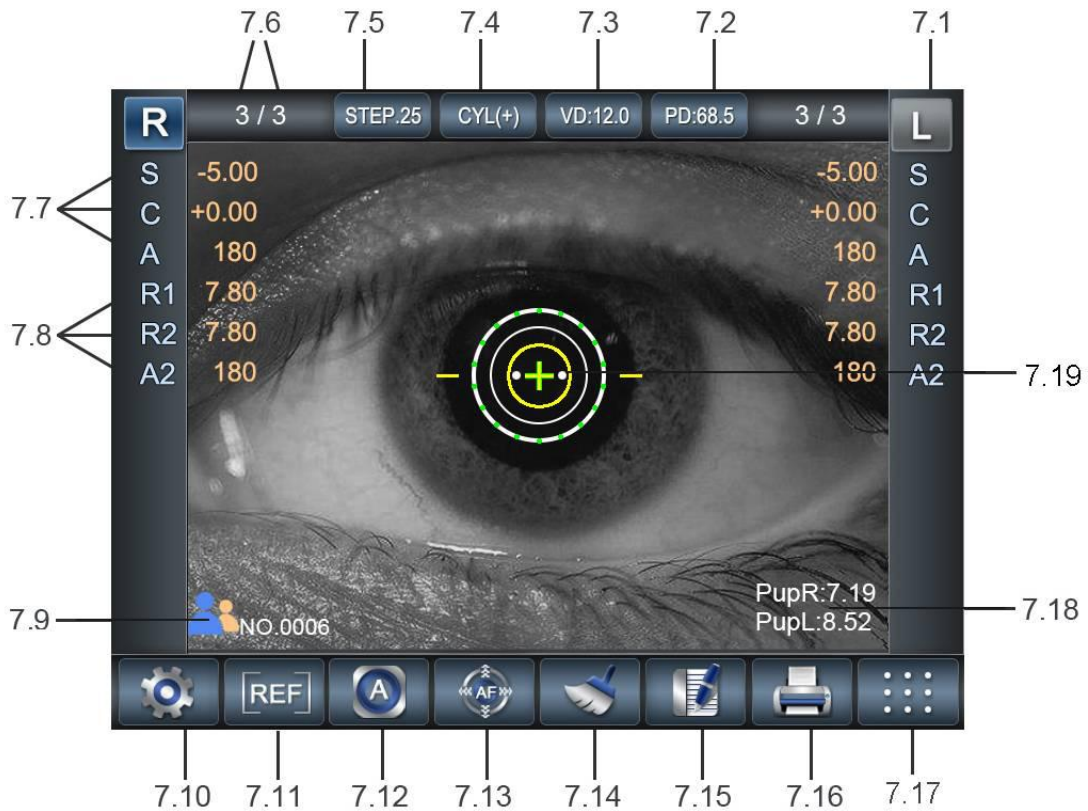

(Figure 7)

- 7.1 L/R Sign: Flashing sign indicates the current measured eye
- 7.2 Pupil distance
- 7.3 VD selection (shortcut key)
- 7.4 Astigmatism symbol selection (shortcut key)
- 7.5 Step selection (shortcut key)
- 7.6 The number of power/corneal parameters measured
- 7.7 Power display
- 7.8 Corneal value display
- 7.9 Adult//child mode selection
- 7.10 Menu set
- 7.11 Measurement mode selection
- 7.12 Auto/manual measurement selection (partial model)
- 7.13 Auto/manual tracking and focusing selection (partial model)
- 7.14 Data clear key
- 7.15 Data record check
- 7.16 Printer
- 7.17 Lattice display
- 7.18 Left/right eye pupil diameter
- 7.19 Pupil alignment target

8. MENU (Figure 8)



8.1 Measurement Mode Selection


Touch this key to pop up three measurement mode menu (as shown in Figure 8), the user can choose to touch any measurement mode menu under need (KER for Keratometry mode, R/K for RefracKeratometry mode, REF for Refractometry mode)



(Figure 8)



8.2 Auto/Manual Measurement Selection (partial model)

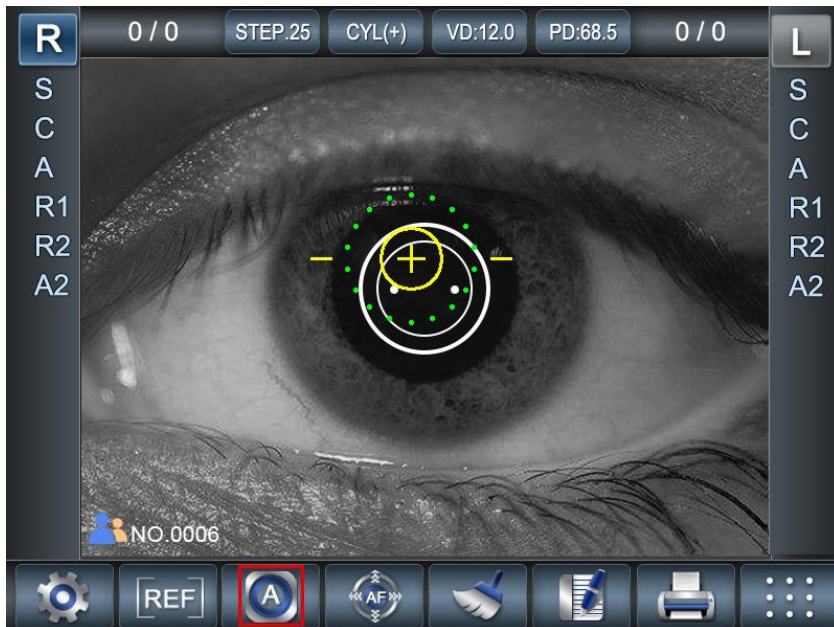Touch ( or ) key (as shown in Figure 9) to select auto measurement mode (A) or manual measurement mode (M)



(Figure 9)

8.3 Auto/Manual Tracking and Focusing Selection (partial model)


Touch( or ) key (as shown in Figure 10) to select auto tracking and focusing


mode (AF) or manual tracking and focusing mode (MF)



(Figure 10)


8.4 Data Clear Key (Figure 11)

Touch  key to clear the measurement data

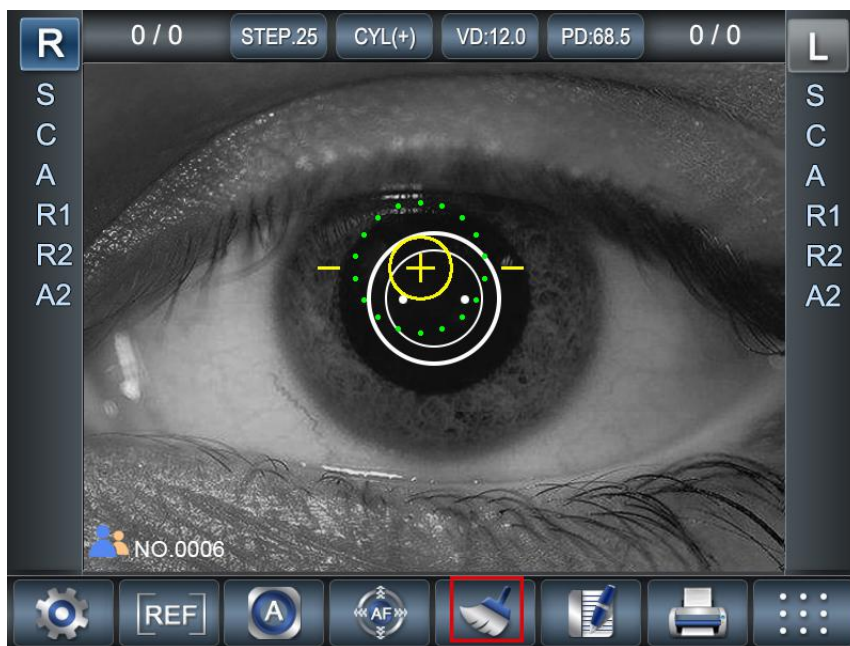

(Figure 11)

8.5 Data Record Check(Figure 12)


Touch ( key to check the measurement data (directly print out the data, the measurement data won't be recorded)




(Figure 12)

Left/right eyes data can be recorded max.10 items separately. Touch REF to display the recorded refractometry data only, touch KER to display the recorded keratometry data only, touch CLEAR key to clear the recorded data, touch RETURN key to return to the measurement interface.


8.6 Printer Set

Touch () key to print current data (Figure 13)

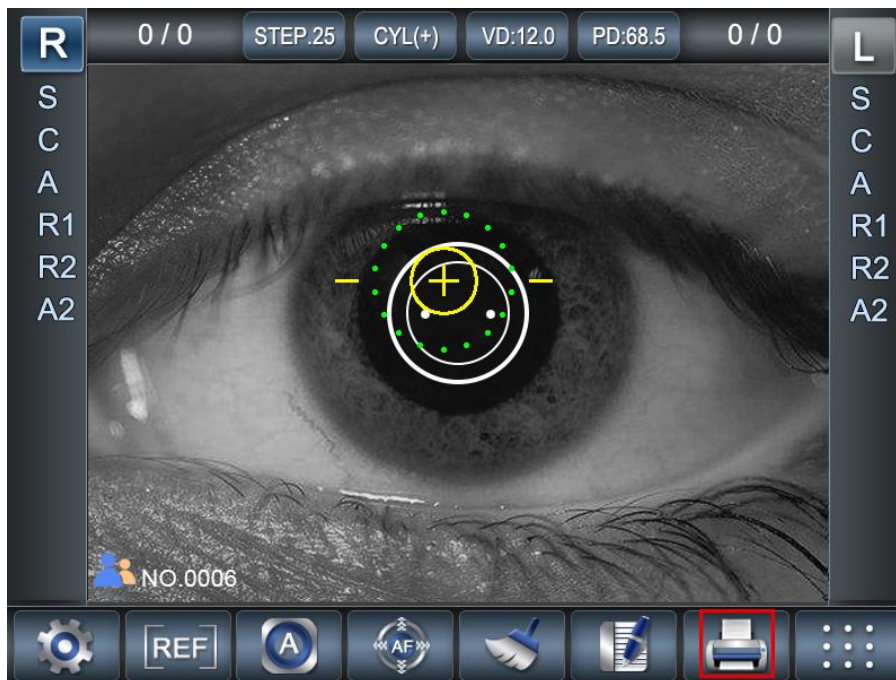


(Figure 13)

8.7 Menu Set

Touch () key to enter the subsidiary menu setting (Current selection for blue)

8.7.1 Refractometry parameters setting (Figure 14)

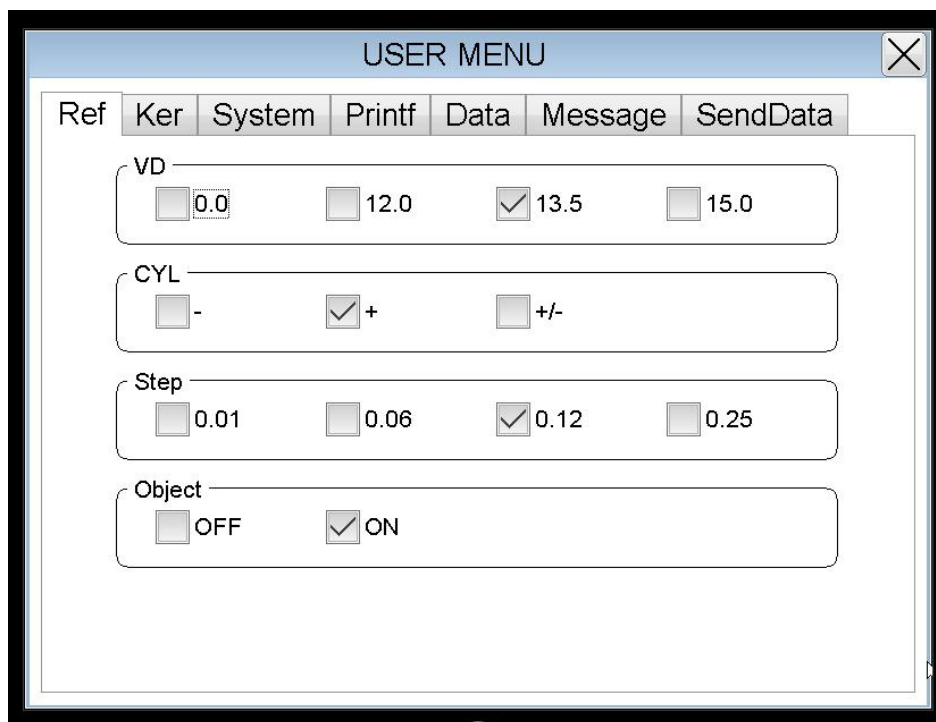

(Figure 14)


VD: Distance between corneal and back top focus of lens, 0.0mm (contact lens), 12.0mm (Asian), 13.5mm (Middle East), 15.0mm (European)

CYL: Astigmatism symbol selection, -, +, ±(Mix)

STEP: Measurement data precision selection

FOGG: Visual guide target atomization function switch (position of guiding target atomization)

8.7.2 Keratometry parameters setting (Figure 15)

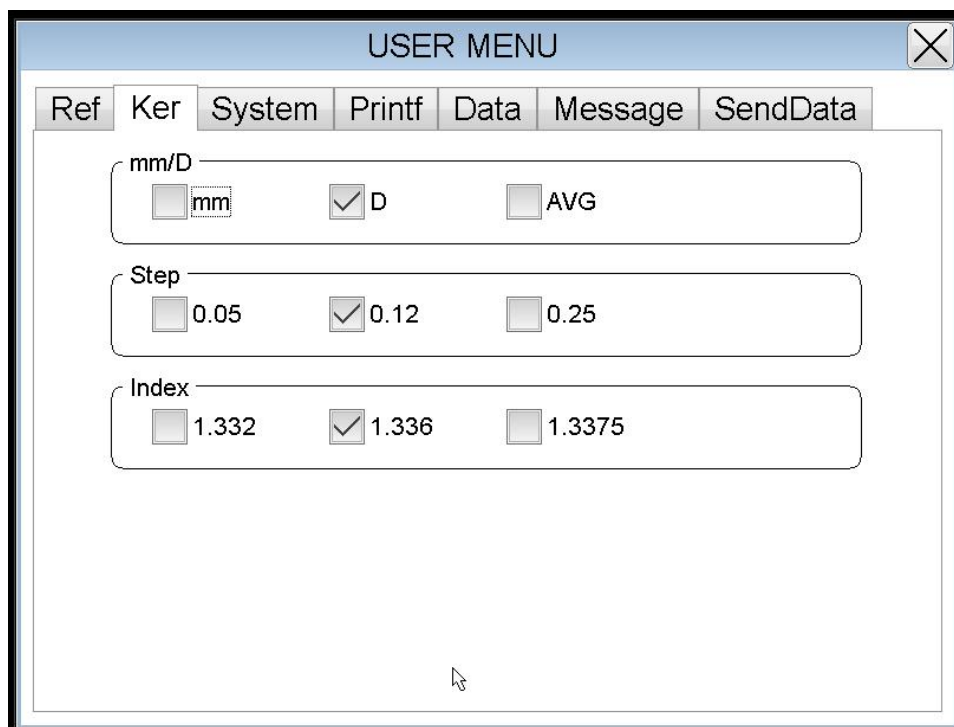


(Figure 15)

MODE: Keratometry radius measurement (mm), keratometry power measurement (m^{-1}) and average value display (AVG) optional

STEP: Keratometry power precision display

REFRACTIVE INDEX: Factory defaults to 1.3375

8.7.3 Mode setting (Figure 16)

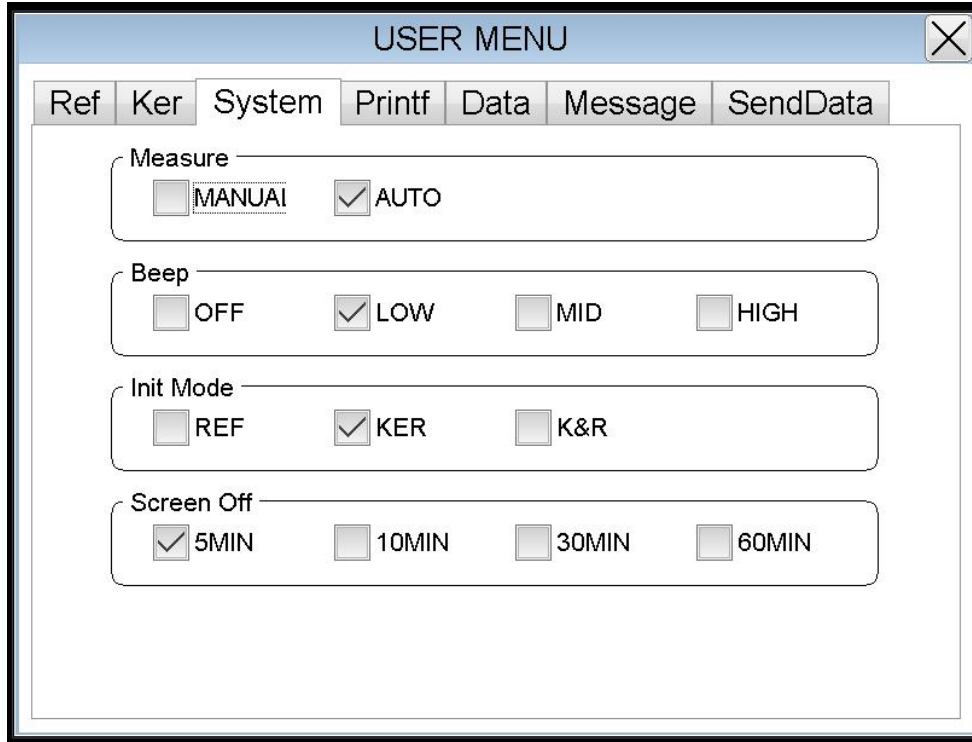

(Figure 16)

MODE: Manual measurement mode and auto measurement mode optional (Auto measurement icon for grey said this model without this feature)


BEEP: Sound prompt when operating. If set off, operation will keep silent

INT-M: Measurement mode selection (same as the main interface function), default startup mode for each starting

STAND BY: Instrument standby time setting (5 minutes, 10 minutes, 30 minutes and 60 minutes optional) (touch any key to wake up)

BRIGHTNESS: LCD screen brightness setting

8.7.4 Printing setting and printing paper replacement (Figure 17)

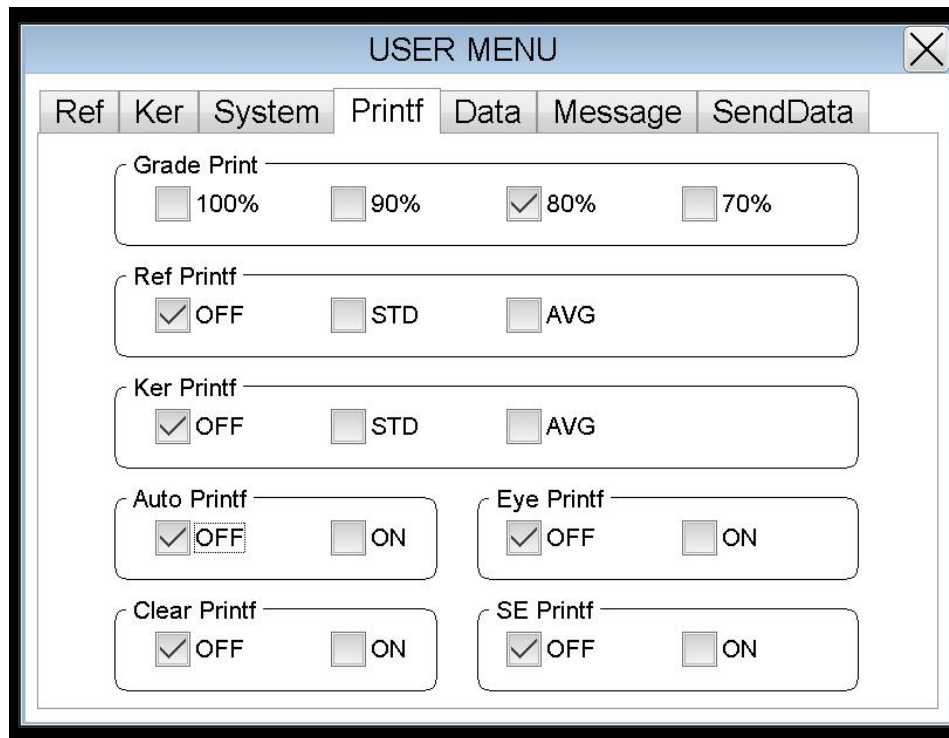

(Figure 17)

PRINT CONCENTRATION :Set the appropriate print concentration according to different thermal printing paper

REFRACTOMETRY: When OFF selected, the refractometry power won't be printed out.

When STD selected, all refractometry power will be printed out. When AVG selected, only print the average value of the refractometry power.

KERATOMETRY: When OFF selected, the keratometry power won't be printed out. When STD selected, all keratometry power will be printed out. When AVG selected, only print the average value of the keratometry power.

AUTO: When ON selected, the measurement results will be printed out automatically after the both eyes measurement finished (in this case, the data is cleared automatically)

When OFF selected, press the print key on panel to print out the measurement results.

EYE: When ON or OFF selected, the refractometry state diagram will be printed or not.

CLEAR PRINT DATA:Automatically erase data after printing


SE: When ON or OFF selected, SE data (the approximate value of cylinder power converted into sphere power) will be printed or not.

PUPIL: When ON or OFF selected, the pupil diameter will be printed or not.

BC: When ON or OFF selected, BC(base curve of corneal contact lens) will be printed or not.

How to install the printing paper (Figure 18)

- 1.) Pull outwards the printer cover, open the cover
- 2.) Place the new printer paper in the box, keep the paper head upwards
- 3.) Draw the paper outwards and directly ride on the printer wheel
- 4.) Push the cover back and close the cover

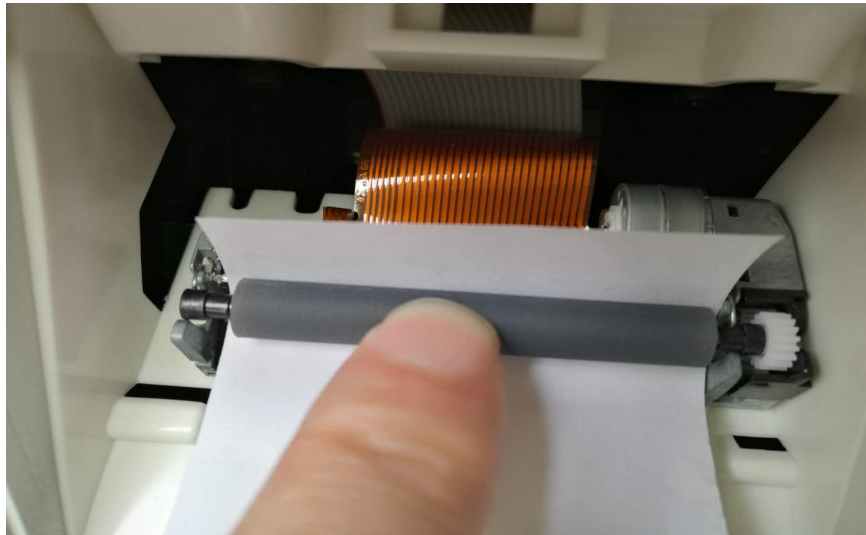

(Figure 18)

8.7.5 Data setting (Figure 19)

USER MENU ✕

Ref
Ker
System
Printf
Data
Message
SendData

Year Month Day

2018 ▼ 1 ▼ 1 ▼

Hour Minute Second

1 ▼ 1 ▼ 1 ▼

No.

0 ▼

Date Format

YMD
 MDY
 DMY

Auto Count

OFF
 ON

(Figure 19)

DISP: Date, month and year display mode

DATE: Edit or modify the exact time of date and month and year

TIME: Edit or modify the exact time of second and minute and hour

COUNT: When ON or OFF selected, recording the number of patients in main interface will be refreshed or not


No.: Patient number setting, patient measuring number setting

Touch DATE, TIME and NUMBER options, enter the sub menu as shown below, select the appropriate number.

8.7.6 Printing message setting (Figure 20)

MSG1 for company name or product model number setting

MSG2 for company address or brand name setting. Users can edit this information freely according to the exact requires. After setting, press ENTER key to preserve and quit.

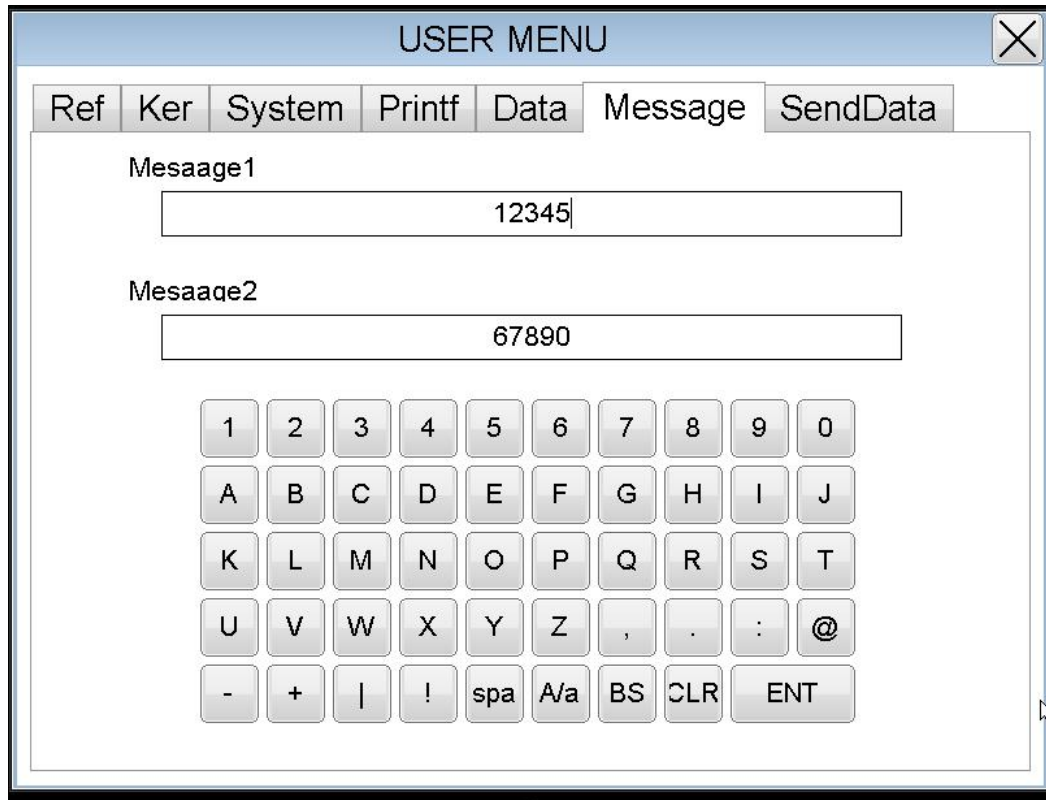

(Figure 20)

ENTER key for confirming and preserving


A/a for capital/small letter conversion

BS for deleting single letter

SPA key for space bar

CRL key for clearing all letters


8.7.7 Data transfer setting (Figure21)

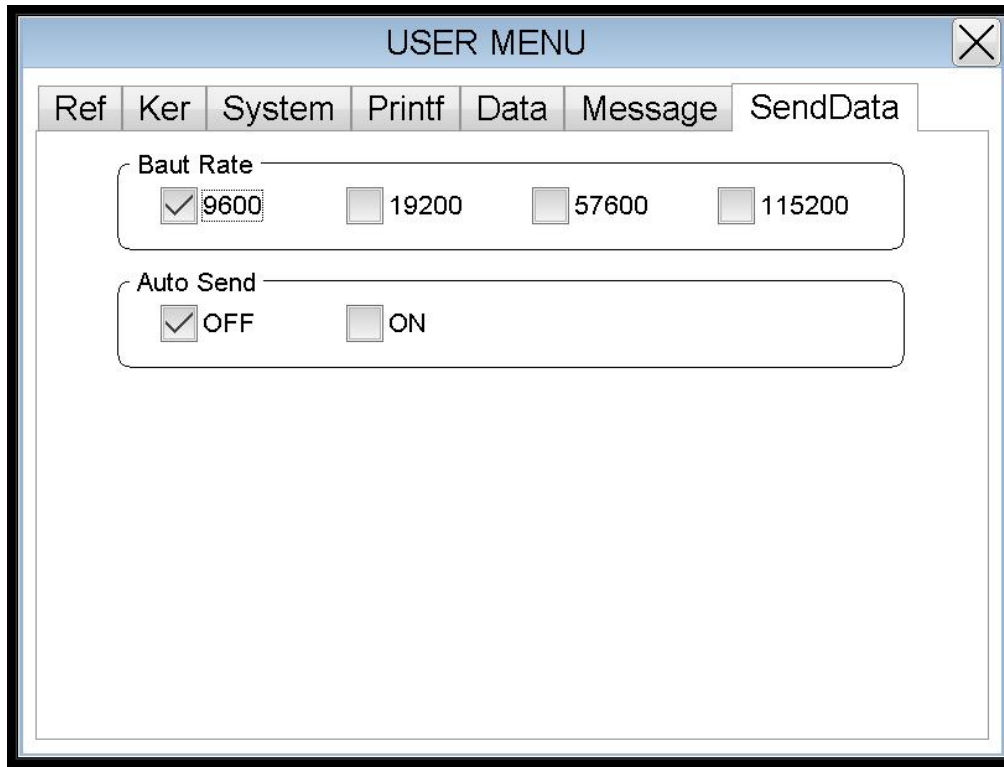

(Figure21)

According to the requires of the connected devices, customers choose the corresponding baud rate, and open the auto option, the measurement data will be automatically transferred to the connected devices, meanwhile the refractometer data will be automatically cleared.

8.8 Array display

touch () key to display the patient's fundus array distribution (indirect evaluation of fundus imaging quality). (Figure.22)



(Figure.22)

8.9 Shortcut Key

8.9.1 Step set: successively touch STEP key to quickly switch 0.01、0.06、0.12、0.25 8.9.2 CYL

axis set: successively touch CYL key to quickly switch —、+、±

8.9.3 VD set: successively touch VD key to quickly switch 0、12、13.5、15

9. MEASUREMENT

Suitable crowd and contraindication

Target patients for adults and children, and crowd of eye power range ($-20 \text{ m}^{-1} \sim +20 \text{ m}^{-1}$). This product is not suitable for newborn eye measurement.

9.1 Preparations before Measurement

9.1.1 Place the device on the specified instrument table, loose the stage fixing lever and keep the device in free sliding state, adjust the four rubber feet to keep the device in horizontal.

9.1.2 Fix and install the specified chinrest paper and printing paper separately

9.1.3 Connect the spare power line to the instrument socket tightly (ensure the local voltage fit to the instrument specification)

9.1.4 Turn on the left side power switch (green indicator light show right in electricity connection), the instrument goes into self-check procedures. After self-check over, it automatically switches to main interface for measurement.

9.2 Notes for Operator and Patient

9.2.1 Adjust the chair height and screen angle in right position

9.2.2 Ensure the patient in comfortable and relaxed posture before measurement

9.2.3 By adjusting the instrument tabletop, keep the instrument height same to the patient natural sitting posture

9.2.4 Settle patient's chin touch the chinrest front and forehead touch the rubber forehead rest in level (keep face parallel with the measurement window)

9.2.5 By observing the patient eyes position and height adjustment mark, press the chinrest up/down key on panel to adjust the patient's eyes same height to the measurement window


9.2.6 By the operation lever, move the sliding body left and right to move the patient's eyes in the measurement range (if the distance of two sides asymmetrical, adjust it by fixing the patient's head deviation)

9.3 Measurement

The measurement alignment method of this device for pupil and center measurement cross target in coincidence

9.3.1 Normal Measurement Mode

Holding the operation lever, quickly shift the sliding body to left side, keeping the measurement window roughly aligning with the patient's right eye socket (Figure.23)

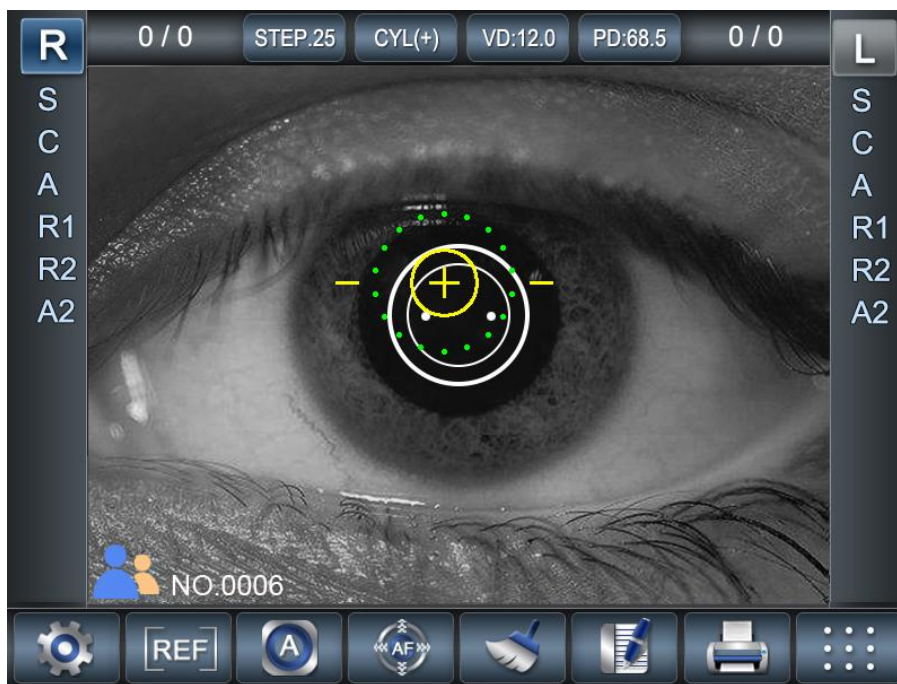

(Figure.23)

Observing the patient's eye location on screen, rotate the operation lever (up and down adjustment), meanwhile swing the operation lever left and right, till the yellow cross-ring target aligning at the patient's corneal vertex, then shift the operation lever front and back, till the patient's eye is clearly focused in the center measurement socket (the accuracy of focusing can be confirmed by observing whether level between the two points of split focusing and cross-ring target) (Figure.24,25)




(Figure.24 away from patient's eye)



(Figure.25 close to patient's eye)

Prompt the patient to open eyes wide (eyelid and eyelash covering eyeball will affect the measurement accuracy), both eyes look right ahead.

Slightly adjust the operation lever, till the two points of slit focusing level with the cross-ring target, and yellow cross measurement target becomes thick and green, press the measurement button, when the measuring light flashing (the screen refreshed in black in moment), it shows the measurement over (the patient no need to see clearly the object-image during measuring, the measurement result same accuracy). The measurement result will be displayed on screen. (Figure.26)。

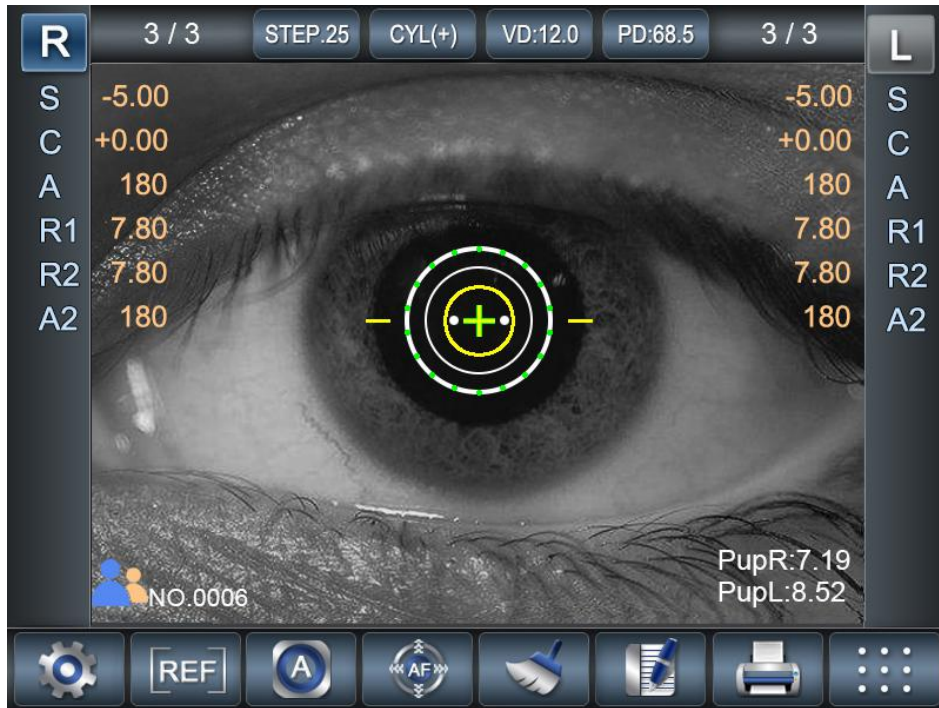


(Figure.26)

Shift the sliding body to right side, repeat the above steps, measure the patient's left eye.

Both eyes measurement over, pupil distance will be displayed automatically on the corresponding position. Choose whether or not to print the measurement results according to settings (auto printing or data output transmission over, the data on screen will be automatically cleared).

9.3.2 Child measurement mode

To measure children or the patients with pupil fibrillation, select child mode (Touch 7.9 key, the right small humanoid icon becomes green). (Figure 27)

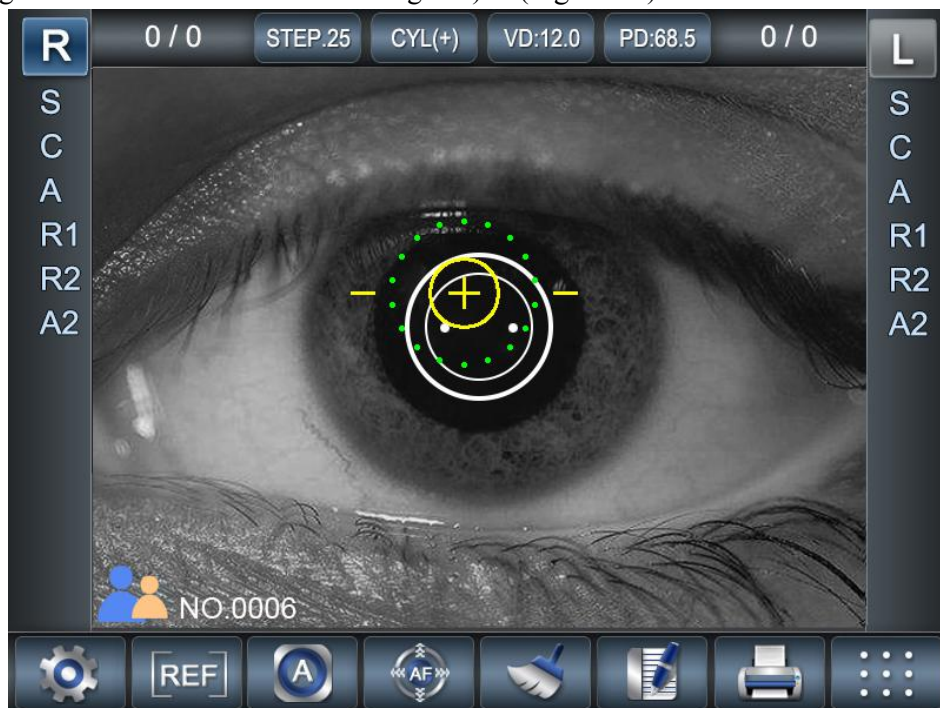

(Figure 27)

9.3.3 Measurement error prompt

During measurement, if the patients found having eyelid ptosis, eyelash disturbance, cataract, microcoria, keratopathy, corneal vertex and pupil center noncoincidence, the error prompt will appear on screen when the instrument can't measure normally, please select the manual measurement mode or force measurement mode (long press the measurement button). (Figure 28 29).



(Figure 28)



(Figure 29)

10. COMMON TROUBLE SHOOTING

10.1 Power indicator light not work

Check and confirm whether the local power fits to the instrument, whether the power plug loose, or whether the fuse damaged (in case this happen, please replace the same specified fuse)

10.2 Chinrest not lift

Check whether the chinrest lift to limit position

10.3 Printer can't work regularly

Check whether the printing paper is finished (in case this happen, red indicator light on panel will flash). Or whether the print setting is correct, and whether there is the measurement data (no data, not print)

10.4 Sliding body not flexible

Check whether the stage fixing lever placed at right position, or whether other sundries go into the slide slot

10.5 Press measurement button, but no data appear

Check whether the patient pupil smaller than 2mm, whether the eye position seriously incorrect, whether the cross measurement target aligns with the patient pupil (the target becomes thick and green), or whether the patient eyeground seriously diseased

10.6 Measurement light not work

When the measurement over, the measurement light will automatically turn off. Sway the sliding body, it will turn on automatically

10.7 If other problems appear, please contact the local agent or original manufacturer

- If the fault phenomenon listed in the common trouble shooting can't be resolved, please contact the original manufacturer or local agent to repair.

- Please provide us with the following information:

Instrument name and model number

Instrument serial number

Fault phenomenon (detailed as possible)

(1.) Accessory maintenance limitation

Providing maintenance accessories to maintain the instrument functions during the instrument lifetime

(2.) Processing of instrument

- To be disposed carelessly of the instrument and accessories will pollute the environment

- Please contact the professional waste disposal company or local dealer before disposing this instrument

11. PACKAGING, TRANSPORTATION, STORAGE

Storage condition between -25°C and +40°C, transportation condition between -40°C and 70°C, relative humidity between 30% and 75%, air pressure between 86kpa and 106kpa

- (1.) Packing list, certificate and manual are included in the packing box
- (2.) The product packaging is not allowed to be shipped with flammable, explosive, corrosive products. Loading should be neat, stable and firm, super high and overweight is not allowed. In transit, rain and snow prevention, anti sun, anti impact, drop prevention should be noted carefully.
- (3.) The product packaging should be stored in a room temperature, dry and well ventilated warehouse, and can't be stored with chemical agents, acid and alkali substances, and other harmful substances.

12. ENVIRONMENTAL PROTECTION

The instruments that have been scrapped, should be strictly disposed in accordance with the requirements of local laws and regulations

13. ENCLOSURE ACCESSORIES

Number	Specification	Quantity
1	User's Manual	1
2	Dustproof Cover	1
3	Lens Dustproof Piston	1
4	Model Eye	1
5	Power Line	1
6	Chinrest Paper	1
7	Cleaning Cloth	1
8	Printing Paper	1
9	Fuse	2
10	Chinrest Pin	2