

Medical Diode Laser System

Reliable, portable & affordable

Built-in battery has been out of date, due to expensive and time-consuming after-sale service. Both power bank and batteries pack are adopted as external power supply


Superpulsed up to 50,000Hz with standard pulse waveform

Wireless foot switch (wired foot switch is optional)

“Hold” & “No-hold” options to activate laser emission

Enjoy music at work

Main applications with various handpieces & tips

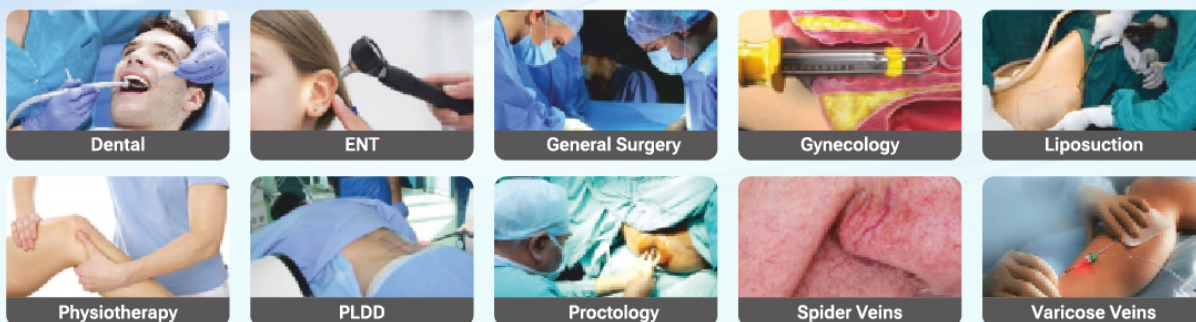

Technical Specifications

Laser Type	Diode Laser		Delivery System	400µm, 600µm, 800µm fibers
Model	Simple 15A	Simple 15B	Aiming Beam	635nm@3mW, adjustable density
Max Power at CW	15W	15W	Voltage	AC 110V-220V, 50Hz/60Hz
Wavelength	980nm	1470nm	Display	8" Full Touch Screen
Operation Mode	CW/Pulsed		Activation	Wireless/Wired Foot Switch
Pulse Duration	10µs-1s		Frequency	Up to 50,000Hz
Dimensions(L*W*H)	280mm*130mm*190mm		Net Weight	2.3kgs