


# Ethylene Oxide Sterilizer

- Effective at low temperature sterilization
- Material compatibility
- Excellent penetration with no lumen restrictions
- Non-corrosive
- Relatively low capital equipment cost



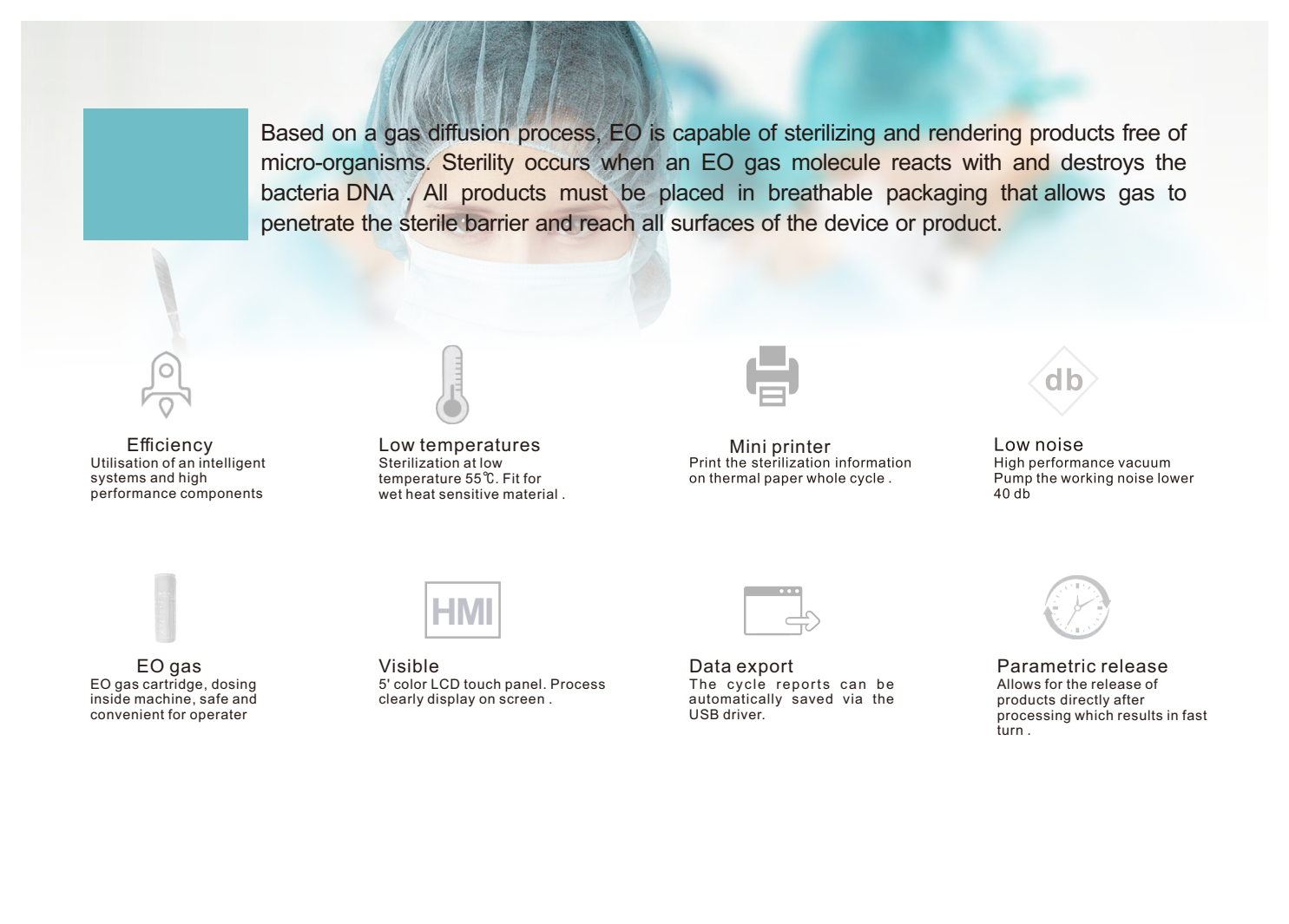

Based on a gas diffusion process, EO is capable of sterilizing and rendering products free of micro-organisms. Sterility occurs when an EO gas molecule reacts with and destroys the bacteria DNA . All products must be placed in breathable packaging that allows gas to penetrate the sterile barrier and reach all surfaces of the device or product.



#### Efficiency

Utilisation of an intelligent systems and high performance components



#### Low temperatures

Sterilization at low temperature 55°C. Fit for wet heat sensitive material .



#### Mini printer

Print the sterilization information on thermal paper whole cycle .



#### Low noise

High performance vacuum Pump the working noise lower 40 db



#### EO gas

EO gas cartridge, dosing inside machine, safe and convenient for operator



#### Visible

5' color LCD touch panel. Process clearly display on screen .






#### Data export

The cycle reports can be automatically saved via the USB driver.



#### Parametric release

Allows for the release of products directly after processing which results in fast turn .

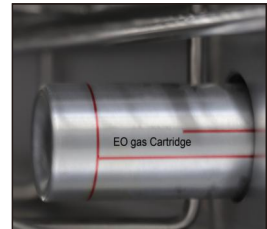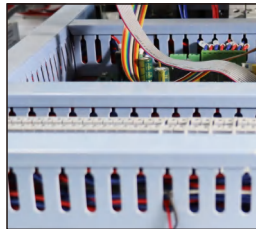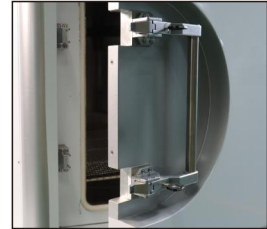


SANQIANG

## Technical parameters

Model	H80	H120	H220	H330
Voltage	AC220V	AC220V	AC220V	AC220V
Power	2.6KVA	3KVA	3.3KVA	3.6KVA
Cell material	304 stainless steel	304 stainless steel	304 stainless steel	304 stainless steel
Chamber size	585X275X422mm	585X375X472mm	705X505X552mm	760X545X712mm
Volume	67L	103L	196L	294L
Machine size	795X615X1140mm	795X615X1140mm	968X800X1240mm	968X800X1340mm
Net weight	148KG	177KG	243KG	295KG
Working pressure	-60Kpa	-60Kpa	-60Kpa	-60Kpa
Cycle time	5-8 hours	5-8 hours	5-8 hours	5-8 hours
Working temperature	30-60°C	30-60°C	30-60°C	30-60°C
Humidity	RH 45-85%	RH 45-85%	RH 45-85%	RH 45-85%
Display screen	5'color LCD	5'color LCD	5'color LCD	5'color LCD
Data export	USB	USB	USB	USB
Mini printer	Yes	Yes	Yes	Yes
Gas cartridge	48 g	72 g	170 g	240 g
Shelves(trays)	1	1	2	2

# EO Sterilization Supplies

---

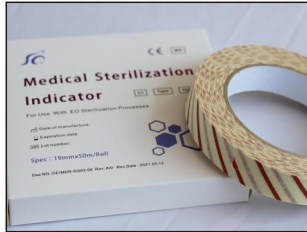

EO indicator tape



EO sterilization pouch



EO indicator strip



EO biological indicator



EO label



Medical sterilization pouch