



Better Detail For Better Quality
Better Detail For Better Quality

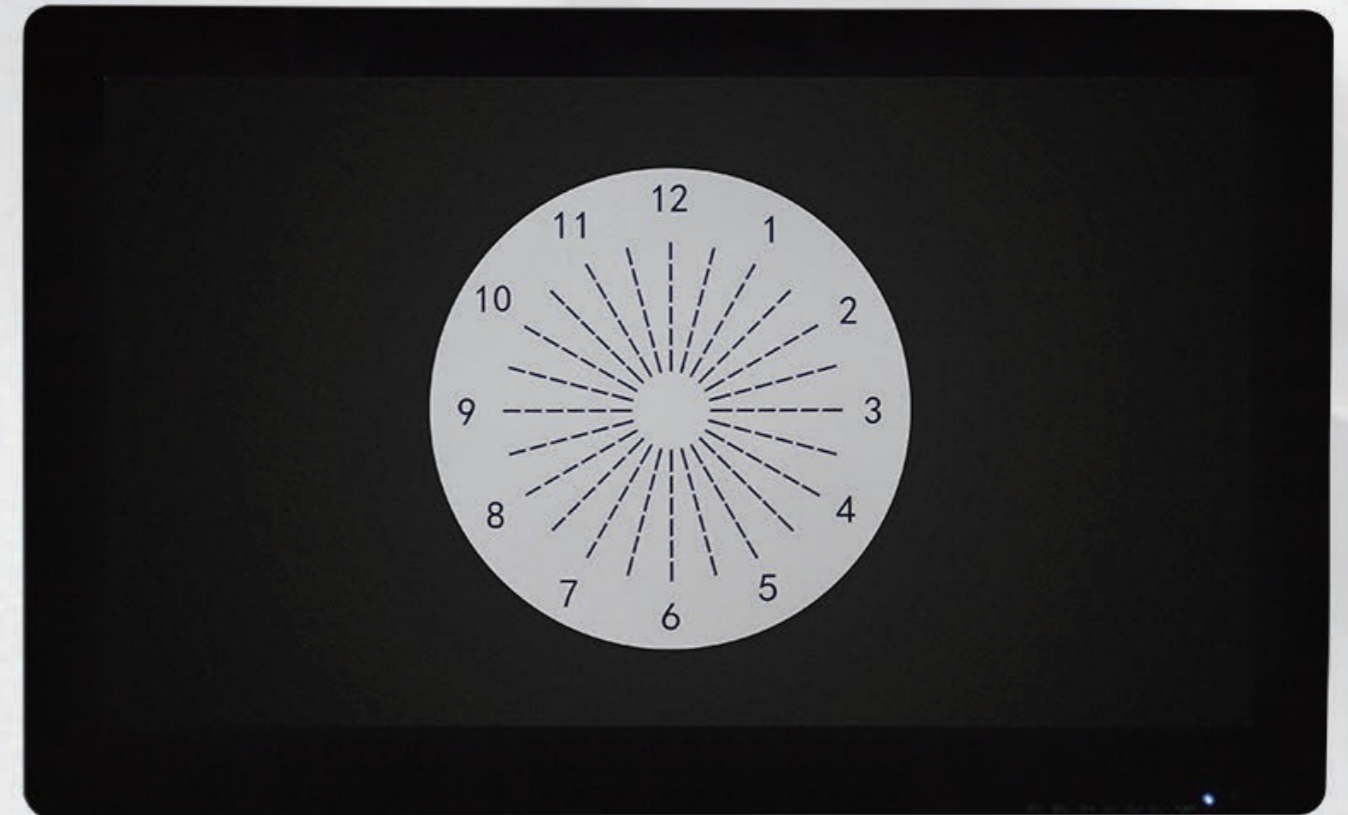

LCD Chart Projector

New Technologies And Powerful Functions

Simple operation interface provides
with accurate operation experience.

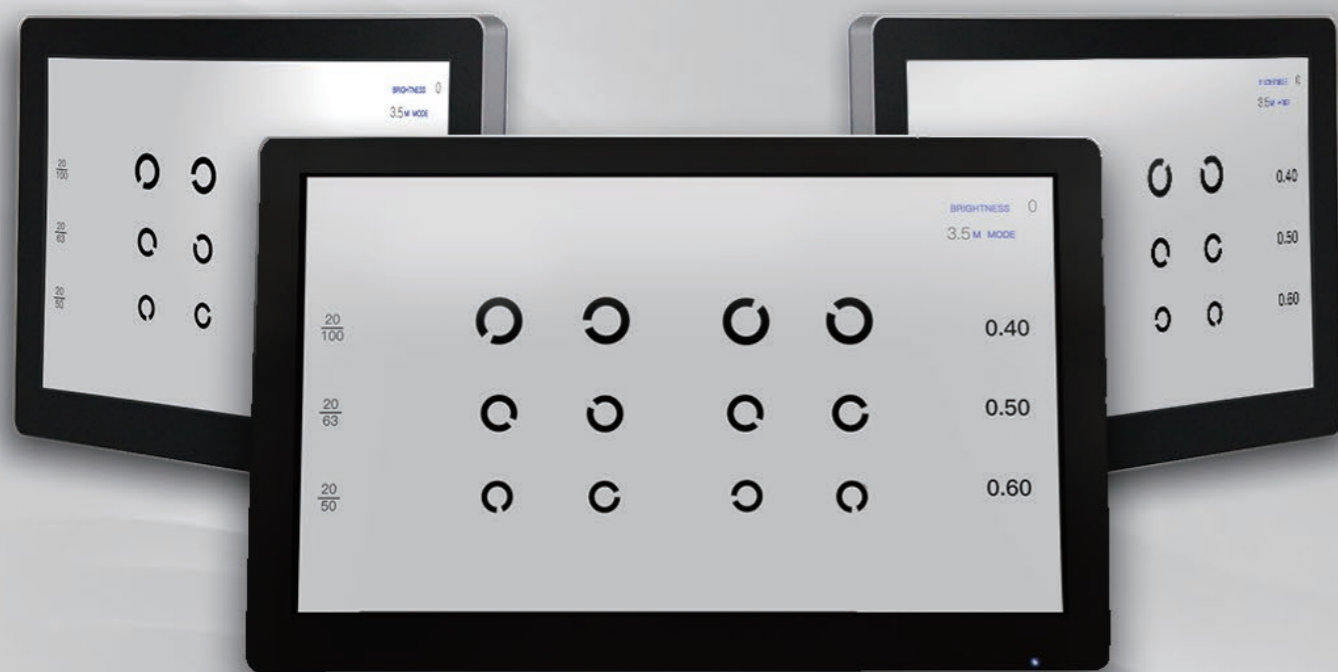

LINK LCD Chart Projector LCP – 150 has multiple extended functions, 23 inch LED backlit screen, polarized charts is provided for binocular test. Various kinds of tests meet the diverse needs of the measurement. Quick test list– Pick ten or less frequently used tests to establish your own quick test list. It can be connected with computer, DVD player or input with multi – media signal through VGA/AV/YPBPR/USB connections as a multi –media player. Excellent programming give full play to the features of CPU.


Connection makes a complete work system

LCP-150 can connect with digital phoropter DVT-1000, which you can control LCD Chart by panel of phoropter




Multi-functional charts meet different testing requestv

E charts, C charts, Letter Charts, Number Charts, Children Charts, Color blindness Charts

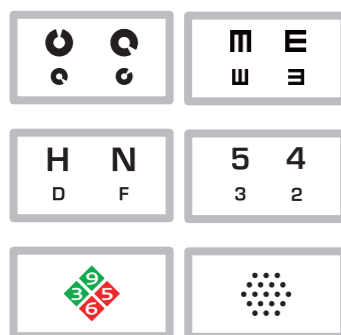


Test distance is adjustable to meet your room space

1–6 meters distance is optional according to room space; it can help to save your space.

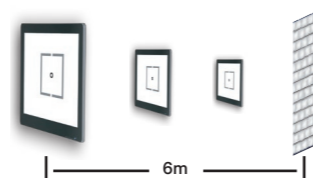

Polarized charts for binocular test

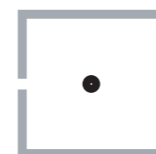


External port on LCP-150 provide for USB flash disk connection

Can run the au-tomatic update procedure
Video playing and image browsing



Easy mounting and operation

Wall and table bracket is provided according to customer's specific request

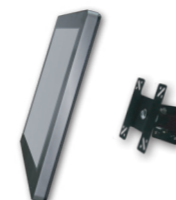

Specification

| | |
|----------------------|--|
| Screen | 23'' LCD (Color) High resolution Monitor |
| | Resolution: 1920*1080 pixels |
| | Dot pitch: 0.2652mm |
| | Screen ratio: 16:9 |
| | Visible Size: 513.8*291mm |
| Brightness | 250cd/m ² |
| Chart Versions | E,C, Letter, Number, Children |
| Working distance | 1–6 meters |
| Remote controlling | Distance > 8m |
| Power Supply | 100–240V AC, 60/50Hz |
| Dimension | 590(W)*365(H)*60mm(D) |
| Weight | 5.5KG() |
| Standard Accessories | Remote Control |
| | Red–green glasses |
| | U Disk |
| | Hanging Bracket |
| | AC–DC Adapter |
| | Power Cable |