


Electric Operating table

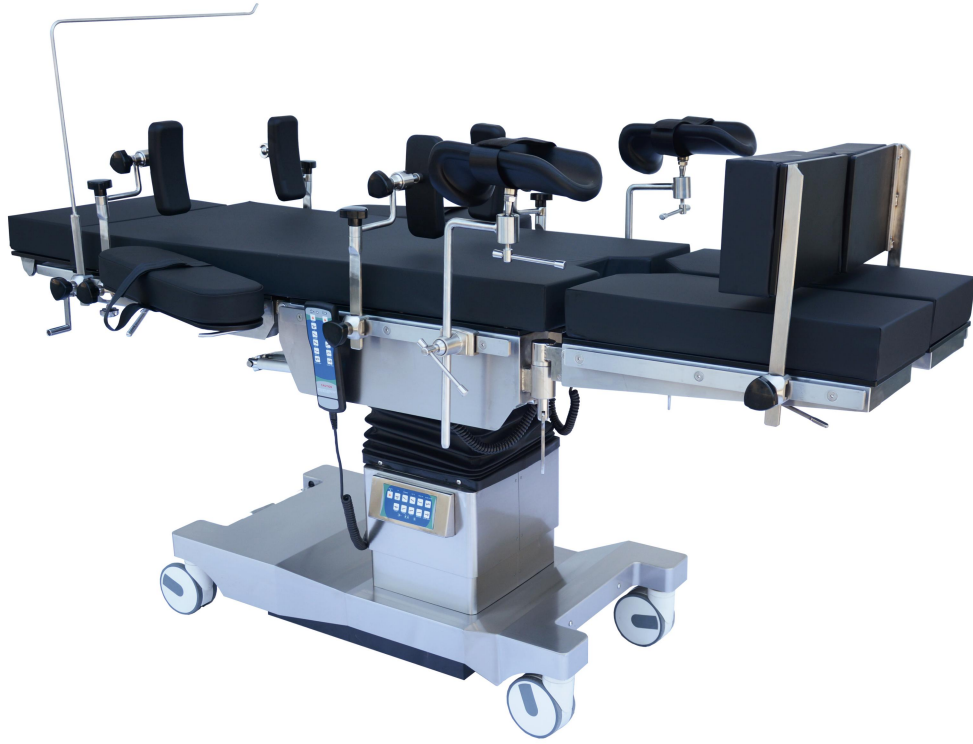

型号 Model	台面长度 Length	台面宽度 Width	台面最低高度 Tabletop lowest height	台面升降行程 Elevation range
	≥2100mm	550mm	≤700mm	300mm

YFD005 can withstand even extreme loads in any position. The structure enables it to be adapted to match the size of the patient and in no time at all it can be equipped for any surgical discipline. A longitudinal shift of up 350mm makes the operating table even more versatile.

Engineered for increasingly complex surgical scenes, YFD005 comes with intelligent automatic control system satisfying various demands of general surgeries. Extra wide tabletop, extreme loading system, long horizontal sliding that can be suitable for both X-ray and C-arm. YFD005 is an ideal electric operating table.

- Longitudinal slide up to 350mm for free access to C arm
- Extreme weight load capacity of 250kgs with highest safety and stability
- Memory pad with antistatic, waterproof design
- Standard override panel control makes the control safer
- Battery power for 50-80 operations

基本参数 Parameter

动作项目 Movements		型号 Model: YFD005
台面前后倾 Tabletop trendelenburg	前倾 Reverse Trendelenburg	$\geq 25^\circ$
	后倾 Trendelenburg	$\geq 25^\circ$
台面左右倾 Lateral Tilt (Left And Right)		$\geq 15^\circ$
头板上下折 Head plate	上折 Up	$\geq 45^\circ$
	下折 Down	$\geq 90^\circ$
	弯曲 Flexible	$\geq 40^\circ$
背板上下折 Back plate	上折 Up	$\geq 80^\circ$
	下折 Down	$\geq 20^\circ$
腿板上下折、外折 Leg plate	上折 Up	$\geq 15^\circ$
	下折 Down	$\geq 90^\circ$
	外折 Outward	$\geq 90^\circ$
腰桥上升高度 Kidney Bridge Elevation		$\geq 120\text{mm}$
台面纵向水平移动 Horizontal Sliding Range		$\geq 350\text{mm}$

标准配置 Disposition

麻醉屏架 Anaesthesia screen	1 piece	支肩架 Shoulder support (cushion)	2 piece	搁手架 Arm support (cushion)	2 piece
支身架 Body support (cushion)	2 piece	托腿架 Leg support (cushion)	2 piece	足板 Foot plate (cushion)	2 piece
长滑块 Long Fixing Clamp	2 piece	滑块 Fixing Clamp	7 piece	手柄 Controller	1 piece
手控器 Hand-held operator	1 piece	床垫 Mattress	1 piece	电源线 Power line	1 piece

包装数据 Packaging and delivery

包装 Packaging	木箱 Wooden case
尺寸 Size	140*74*103cm
毛/净重 G./N.W.	290/240kg