




INTELLIGENT MICROPLATE WASHER

Features

1. Microcomputer control, 5.1-inch LCD big screen; All wash information are displayed in a screen.
2. Humanized menu, Humanized-machine interaction mode, convenient operation and no need of training.
3. Apply to multiple ELISA Plate, such as plain, U, V and round type.
4. By replacing the washing head (8/12) to portrait or landscape washing
5. Can randomly edit and store 100 group programs, and display all the parameters on the same screen.
6. Have the functions of vibrating plate and soaking, the time is adjustable.
7. Have the functions of washing the base and three-point inhaling the liquid. Ensure the washer residues to a minimum that the cleaning effect is better.
8. There is a design of two plate incubation place with inner-circulation wind shower.(optional)
9. Incubation time can be adjusted between 0~90minutes.(optional)
10. Alarm automatically when lacking of lotions or the waste liquid full.
11. The full waste liquid to be discharged automatically. (optional)
12. The design of the transparent shield in the area of washing can screen the noise, prevent the spilling and spatter, and guarantee the safety.
13. Four patent technology.

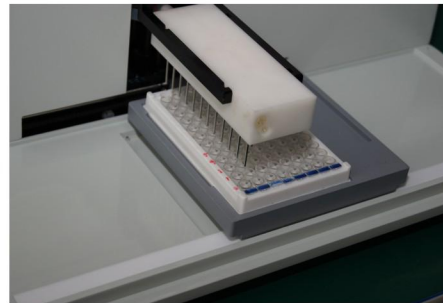

Specification:

1. The head of washer: 8 /12pin
2. Residual Volume: <0.8ul/well
3. Dispense precision: CV<1.5%
4. Dispense Volume: 50–5000ul, adjustable
5. Washing times: 0~900 times, adjustable
6. Soaking time: 0~9000S
7. Vibration time: 0~9000s
8. Memory: up to 100 programs
9. Incubation: 2 plates wind circulation incubation (optional)
10. Temperature: 25–45°C ± 5°C(optional)
11. Incubation time: 0–90mins adjustable (optional)
12. Wash bottle: three 2.5L
13. Interface: Serial port RS–232
14. Waste from the row: Waste liquid full automatic discharge (optional)

China Care Medical Equipment Co.,Ltd
 Address: Room B101, Zhongbang Business Building, PuXing
 Road No. 87, Tianhe, Guangzhou,
 Chinaweb: www.chinacaremedical.com