



Anal examination system
(Optical Endoscope Camera System)
CCD-200

Features


- > High-strength corrosion-resistant aluminum alloy cart. Ultra-quiet self-locking casters with lifting function. 2-axis rotating display stand.
- > Imported footswitch: IPX8 rated.
- > Digital HD handle camera: built-in light source (turns on when in use), automatic white balance. 2 wireless mouse and keyboard.



introduction independently developed endoscope camera system software.

> Basic Functions

- 2 Real-time dynamic display of images, image acquisition, and storage functions.
- 2 Freeze frame.
- > User Profiles: User profiles can be created, saved, selected, and deleted through the user menu. > Print Function.

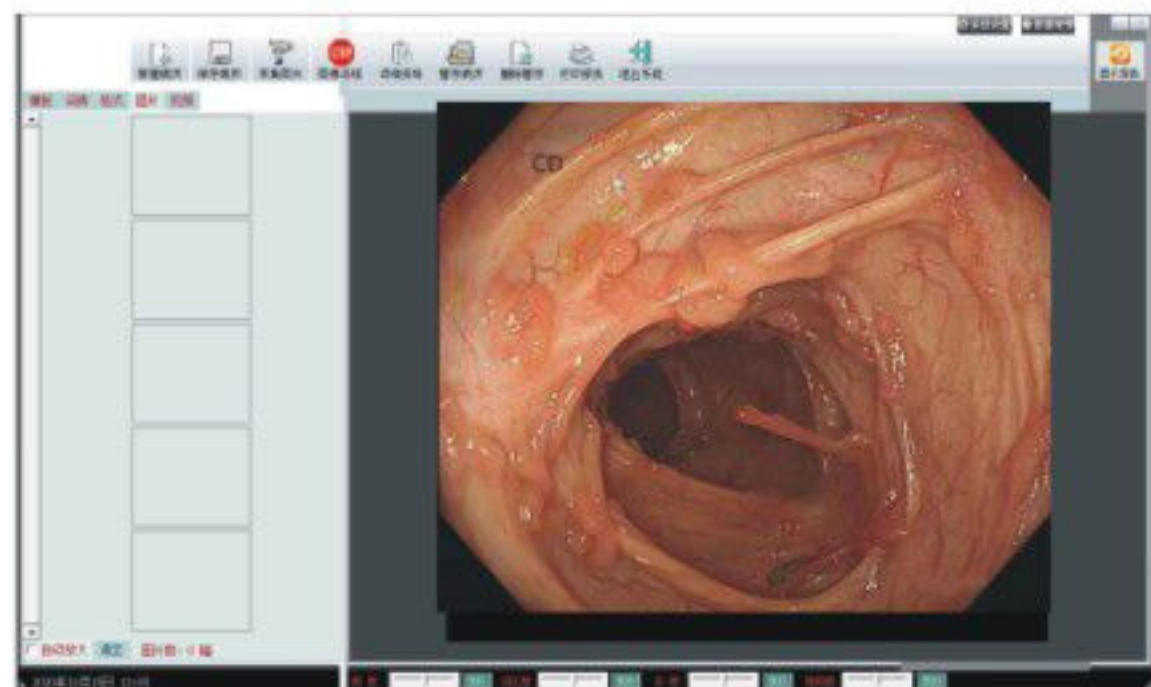

> Image Processing Function 2

- | | | |
|-----------------------|-------------------|----------------------------|
| Brightness Adjustment | 2 Mirror | 2 Image Clarity Adjustment |
| Contrast Adjustment | Image Enhancement | Image Measurement |



> Software Security Protection

The software dog protects patient records, and the user can only obtain permissions for system activation, operation, data file management, etc., when the software dog is in use.

Digital HD Handle Camera

- Horizontal resolution greater than 800 lines, vertical resolution greater than 800 lines, high resolution.
- > USB digital signal output, image clarity.
- > Compact size, portable, easy to operate.



Can be equipped with disposable anal speculum, available in large, medium, small, and pediatric sizes, safe and reliable.

Endoscope camera and light source (optional)

► Digital HD portable endoscope camera

- A USB data 3.0 output
- Supports one-key white balance
- Can be equipped with disposable rectal speculum and disposable sigmoid colon speculum, allowing for accurate examination and diagnosis of the rectum and sigmoid colon.



China Care Medical Equipment Co.Ltd

Room B101, ZhongBang Business +86 139 2648 0041
 Building, PuXing Road No. 87, Tianhe info@chinacaremedical.com
 Guangzhou, China