

Traditional biological microscope

A pair of sharp eyes
opens the door to the
micro world for you



Excellent Functions



Pressing Light Source
Adjustment Knob



Light Source State
Display



Biological + Stereo



Portable Handle

Basic Functions



Low voltage power
supply ensures the
safety



Optimal design of
gravity center avoids
the falling down




Bidirectional mobile
stage controls accurate
displacement



Optional battery for
using at any time and
anywhere

Multi-function is integrated into one button and a single touch

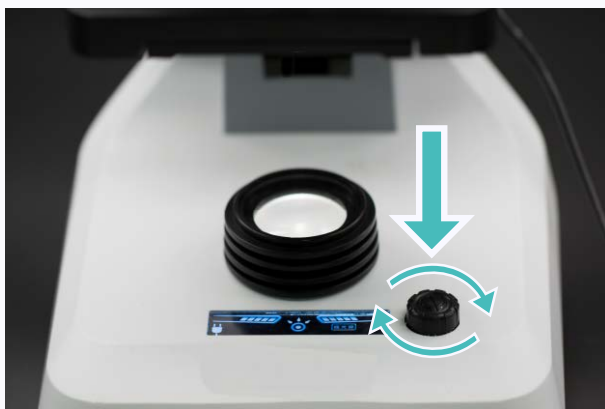

- Long pressing controls on & off
- Short pressing controls light source




Bottom light
source



Side light
source

- Rotary knob controls brightness

Light Source & Electricity State Display

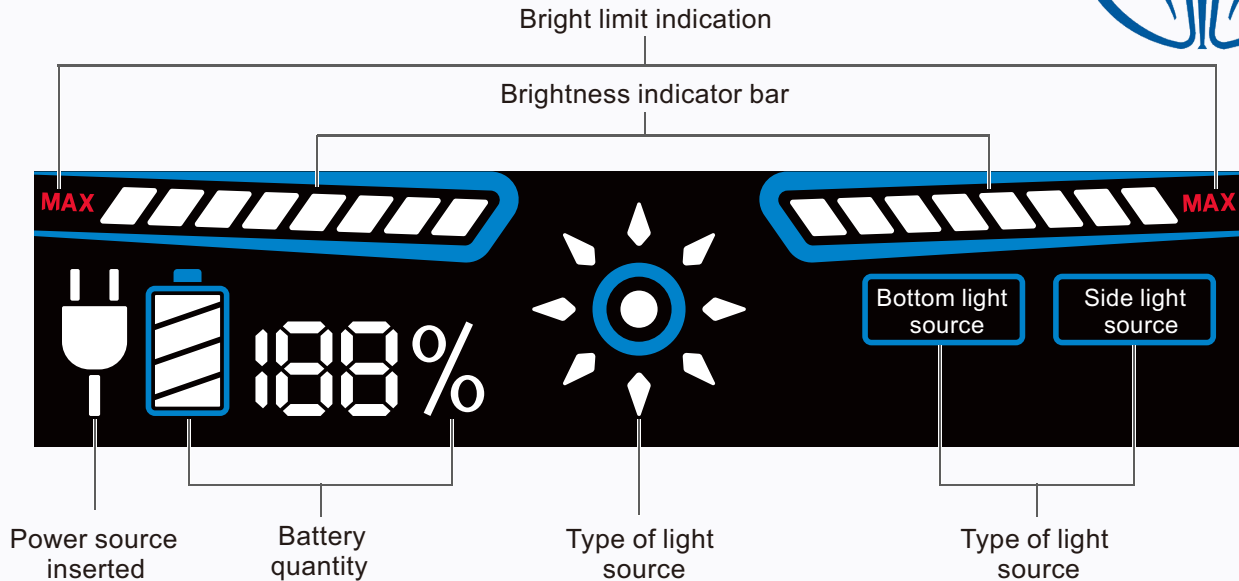

Display on the instrument equipped with a rechargeable battery only

Observation mode is not single one machine dual functions

Besides a bottom light source, it also adds a lateral reflection light to provide a wide range of physical lighting, which makes the microscopes observe not only biological loading but also physical entity.



Biological



Stereo



High capacity rechargeable battery allows you to observe without interruption.


The instrument can be equipped with a high capacity lithium battery(optional) , and can be used without connecting the power supply, which is convenient for outdoor observation anytime and anywhere.

Microscope in hand between coming and going

The eye-catching handle on the top of the instrument defines the holding district , which insures the microscope is upright when lifted, avoids the falling down of the eyepieces and the sample and effectively prevents the loss of accuracy.



SPECIFICATION

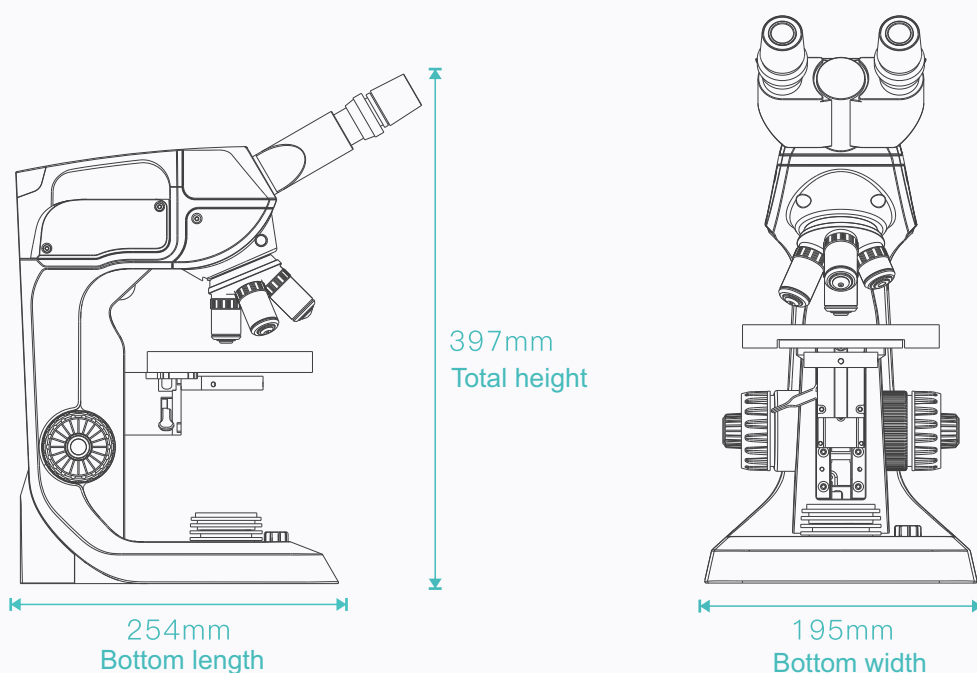

Model	CCL-AC1	CCL-AC2
Type	Monocular	Binocular
Microscope function	scoBiological & Stereoscopic Micrope	
LED lighting	Bottom & Side light source	
Focusing handwheel	Coaxial Coarse & Fine Focus System	
Stage	Double Layers Mechanical Stage	
Condenser	Plug-in Condenser	
Input Voltage	19V	
Eyepiece	WF10×/ 18mm (Optional WF10×/ 20mm)	
Converter	3 - Hole	4 - Hole
Biological Objective	Finite 4X/10X/40X	Infinity 4X/10X/40X/100X
Stereoscopic Objective	Optional	
Optical magnification	40X – 400X	40X – 1000X
Packing List	Power Cord, 19V Adapter, Dust Cover, Micro Slide, Warranty Card,	
Optional Accessories	Aluminum alloy box/Aluminum alloy box with bar/High capacity lithium battery	

China Care Medical Equipment Co., Ltd.
Address: Room B101, Zhongbang Business Building, No. 87 Puxing Road, Tianhe, Guangzhou, China Web: www.chinacaremedical.com